CONTACT INFO

✉ magazzino@evolutamente.it

🌐 www.elioslamp.com

🌐 www.evolutamente.it

USER MANUAL

NOVA BRUSH


SONIC ELECTRIC
TOOTHBRUSH
WITH LED TECHNOLOGY

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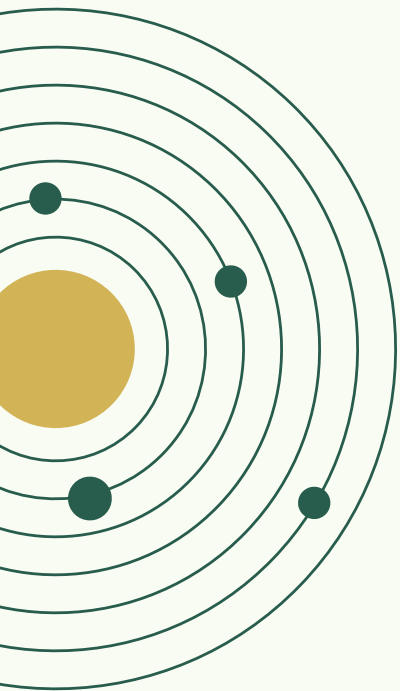USER MANUAL

NOVA BRUSH

SONIC ELECTRIC TOOTHBRUSH
WITH LED TECHNOLOGY

Before using the device, please carefully read the manual to know all the functions to ensure safe use.

1. Electrical Safety: Never use near water, do not touch the power plug with wet hands, and do not charge in the bathroom. If the charger falls into water, unplug immediately without touching it.
2. Do Not Modify: Do not disassemble, modify or repair the product. This may cause fire, electric shock or injury. Contact a service center for repairs.
3. Proper Use Only: This product is designed to clean teeth, gums and tongue only. Discontinue use and consult a dentist if discomfort or pain occurs.
4. Keep Away from Children: This product is not intended for use by children. Keep out of reach of children.
5. Maintenance: Replace the brush head every 3 months. Do not use if the product, cable or plug appears damaged in any way.

TABLE OF CONTENTS

ENGLISH

<i>How it works</i>	1
<i>Technical specifications</i>	3
<i>What's in the box</i>	4
<i>Operating modes</i>	5
<i>Usage Protocols</i>	7
<i>Maintenance and care</i>	8
<i>Warranty and support</i>	9
<i>Problem Solving</i>	10

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1 HOW IT WORKS

EN

LED phototherapy plays a fundamental role in advanced oral hygiene, influencing not only surface cleaning but also deep oral health, bacterial control and overall gum well-being.

The dual interchangeable head system allows you to choose the most suitable LED therapy for your oral needs, combining advanced sonic vibration with targeted phototherapy.

RED LED HEAD

Anti-inflammatory and Regenerative Therapy

Red light triggers a photochemical reaction that:

- Improves blood circulation in the gums
- Strengthens local immune response
- Stimulates collagen production
- Accelerates gum tissue recovery
- Reduces inflammation and swelling

Ideal for: gingivitis, bleeding, sensitivity, post-dental treatment recovery

BLUE LED HEAD

Antibacterial and Whitening

Therapy Blue light acts on:


- Elimination of oral pathogenic bacteria
- Prevention of plaque formation
- Control of microbial proliferation

- Whitening effect on surface stains
- Maintenance of oral hygiene

Ideal for: daily hygiene, cavity prevention, whitening, fresh breath

HIGH FREQUENCY SONIC VIBRATION

With 38,400 movements per minute, sonic technology completes the LED phototherapy action by effectively removing plaque, reaching interdental spaces and stimulating gum microcirculation.

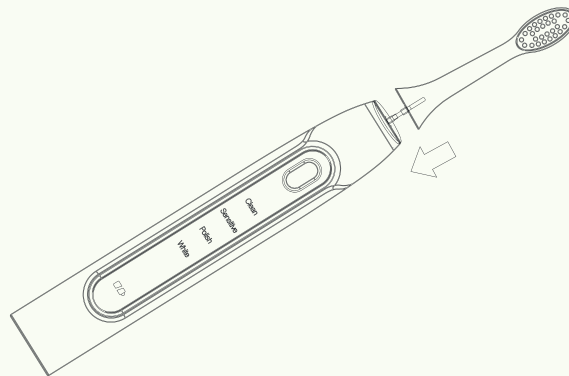2 TECHNICAL SPECIFICATIONS

LED TECHNOLOGY	Red Light + Blue Light
LED PER HEAD	6 LED
VIBRATION FREQUENCY	38,400 rpm
BRUSHING TIME	2 minutes (2x/day)
OPERATING MODES	4 (Clean, Sensitive, Polish, White)
CHARGING VOLTAGE	5V
BATTERY LIFE	Up to 4 weeks (2x/day)
CHARGING PORT	USB

3 WHAT'S IN THE BOX

EN


1. MAIN BODY

2. INTERCHANGEABLE HEADS

3. USB CABLE

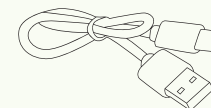4. CHARGING BASE


4 OPERATING MODES

To select a mode, press the button until the desired indicator lights up.

Each mode lasts 2 minutes with pauses every 30 seconds to guide cleaning of the 4 mouth quadrants.

The memory function retains the last mode used.

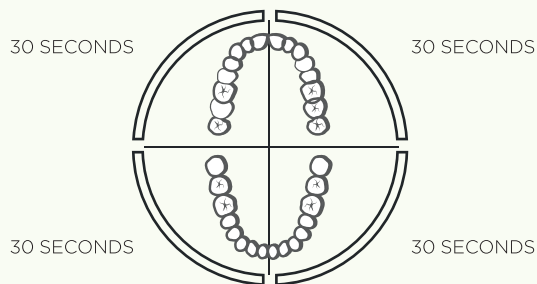BUTTON FUNCTIONS GUIDE

1 CLEAN - Standard Cleaning

Basic mode for effective plaque removal. Perfect for daily morning and evening use.

Recommended combinations:

- BLUE Head: daily hygiene and prevention
- RED Head: after dental treatments or for delicate gums.

2 SENSITIVE - Gentle Mode

Reduced vibration for sensitive teeth and gums. Ideal for those who find sonic vibrations too intense.

Recommended combinations:

- RED Head: acute gingivitis, bleeding, inflammation
- BLUE Head: dental sensitivity, first weeks of use

3 POLISH - Quick Polishing

Rapid oscillating movement to polish tooth surface. Optimizes whitening effect and brightness.

Recommended combinations:

- BLUE Head: extra treatment for incisors and visible teeth.
- Use once daily, preferably in the evening

4 WHITE - Intensive Whitening

Maximum power to remove stubborn stains from coffee, tea, red wine and smoke.

Recommended combinations:

- BLUE Head: 4-week whitening program.
- Use with whitening toothpaste or gel for optimal results.

5 USAGE PROTOCOLS

SYMPTOMS	HEAD	MODE
Daily hygiene	Blue	Clean
Gingivitis	Red	Sensitive
Inflammation	Red	Sensitive
Maintenance	Blue + Red	Clean
Intensive whitening	Blue	White + Polish
Gum sensitivity	Red	Sensitive
Dental sensitivity	Blue	Sensitive

6 MAINTENANCE AND CARE

EN

FIRST USE

1. Fully charge the battery (approximately 4 hours)
2. Choose the appropriate head for the desired protocol
3. Insert the head onto the handle until you hear the click
4. Wet the bristles and apply toothpaste

CORRECT USE

1. Position the bristles at 45° on the gum line
2. Press the button to turn on and select the mode
3. Gently guide the toothbrush without pressing hard
4. Pauses every 30 seconds indicate when to change quadrant
5. After 2 minutes the device turns off automatically

CLEANING AND MAINTENANCE

1. After each use: rinse head and bristles under running water
2. Once a week: clean the LEDs with a slightly damp cotton swab
3. Do not immerse the handle in water
4. Dry the heads after use and store vertically

BATTERY CHARGING

1. Use only the supplied charger (5.7V)
2. Flashing red indicator = charging
3. Solid blue indicator = fully charged
4. DO NOT use during charging

7 WARRANTY AND SUPPORT

Dear Customers,

Thank you for choosing our product. Customer satisfaction is our primary commitment. For your peace of mind, we offer:

- A limited 2-year warranty from the date of purchase
- Warranty coverage starts from the purchase date indicated on your invoice
- Please keep the original invoice for any warranty claims

WARRANTY EXCLUSIONS

The warranty does not cover:

- Damage caused by improper operation or use
- Malfunctions caused by disassembling the device
- Malfunctions or damage caused by accidents or human activities (scratches, water immersion, or unauthorized parts replacement)
- Damage caused by natural disasters

For technical issues, malfunctions or warranty requests, please contact:

Email: magazzino@evolutamente.it

Avoid Rain

Fragile

Up

Disposal

CE Mark

Avoid Heat

8 PROBLEM SOLVING

EN

FAULT	CAUSE	SOLUTION
Device does not turn on	• <i>Battery discharged</i>	1. <i>Hold power button 2 sec 2.</i> 2. <i>Fully charge</i> 3. <i>Turn on after charging</i>
No light emission or weak emission	• <i>Head not inserted properly</i>	1. <i>Turn off device</i> 2. <i>Reinsert until click</i> 3. <i>Try again</i>
Cracks in body or treatment window	• <i>Physical damage (drops, impacts) or internal malfunction</i>	1. <i>Stop use immediately</i> 2. <i>Contact customer service</i> 3. <i>DO NOT attempt to use</i>

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